

ALPS™ cosmetic covers are pre-stretched and are heated to shrink to size. After pulling skin over prosthesis, use a heat gun to start shrinking the skin at the toes and work up to the proximal brim. If needed, a vinyl contact cement can be used or super glue for a permanent bond. Can be used on transtibial patients. For use with suction suspension applications, place a small hole so the cosmetic cover does not inflate.

- Highly conformable
- Two sizes to accommodate most users
- Pre-Stretched, heat shrink to fit
- No transition line at the ankle

PRODUCT NO.

TSK-X-L	EasySkin™ Large
TSK-X-S	EasySkin™ Small

Color Swatches are available.
COLORS: 3,5,7,9,11,12,13

US Patent No. 6,740,124



Made in the USA




ALPS®

Tel: 727.528.8566 Tel: 800.574.5426 www.easyliner.com
info@easyliner.com

EXTRAS

EasySkin™ TSK Cross Section Diagram

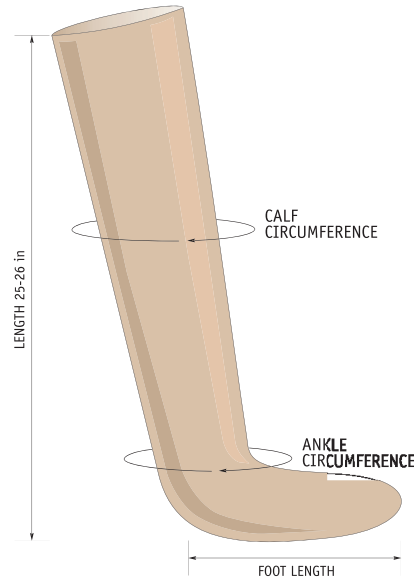

EasySkin™ TSK Sizing Chart
Sizing Specifications

	Foot Length	Ankle	Calf
SMALL			
Minimum	18 cm	17 cm	25 cm
Maximum	23 cm	22 cm	35 cm
LARGE			
Minimum	23 cm	22 cm	35 cm
Maximum	30 cm	29 cm	45 cm

*All ankle and calf measurements are in centimeters.

